


SMARTSCOPE™

ROUNDED TIP LASER FIBER




Medical Device ACCESSORIES

The ForTec Smartscope™ Laser Fiber is compatible with Holmium and Thulium Fiber Lasers. It features a rounded tip which ensures it passes smoothly through the working channel of a fully-flexed ureteroscope. It comes conveniently packaged sterile for single-use.

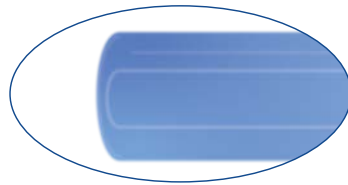

A SmartScope, rounded tip fiber

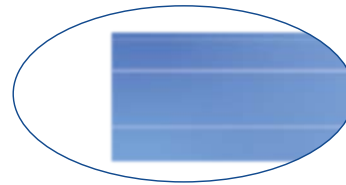

A traditional flat tip fiber

- **Smart Design:** The smooth rounded tip passes effortlessly through the scope's working channel reducing damage and costly scope repairs compared to standard flat tip fibers.
- **Improved Performance:** The fiber tip is hardened resulting in a highly durable tip with more consistent energy transmission.
- **Streamlined Operation:** The fiber is equipped with RFID technology enabling automatic recognition of the diameter and self-adjustment of emission settings.

Fiber Type	Max Input Power	Fiber Core Size	Outer Diameter
ForTec TRUE 200 micron laser rounded tip fiber	25 Watts	200µm	375µm
ForTec 272 micron laser rounded tip fiber	50 Watts	272µm	400µm
ForTec 365 micron laser rounded tip fiber	105 Watts	365µm	550µm
ForTec 550 micron laser rounded tip fiber	105 Watts	550µm	750µm
ForTec 1000 micron laser rounded tip fiber	105 Watts	1000µm	1400µm

[THE FORTEC ADVANTAGE]

- ✓ Service and technical expertise
- ✓ Experienced, professional technician support
- ✓ A stable, professional partner
- ✓ Large inventory of surgical technologies

Call today

800.963.7101