



Case-by-Case
MOBILE SERVICE

neoV Laser

The neoV laser is a remarkably small, lightweight and portable unit that delivers a powerful and safe laser beam for colorectal procedures.

Fistula and Pilonidal Sinus Treatment

The Corona Probe technology distributes the laser energy circumferentially and uniformly, facilitating efficient thermal ablation with minimal risk to the underlying muscle structure.

Coagulation of Hemorrhoids

The laser provides the ability to coagulate and treat hemorrhoids of grades 2 - 4 in a minimally invasive way.

- Proprietary cooling technology
- Energy delivered at 1470 nm
- Specially designed fiber tip with small outer diameter
- Dedicated laser marking



[THE FORTEC ADVANTAGE]

- ✓ Avoid capital and maintenance costs
- ✓ Experienced, professional technician support
- ✓ A stable, professional partner
- ✓ Large inventory of surgical technologies

Call today

800.963.7101

to start ForTec Mobile Service

neoV Laser

Laser Specifications:


Laser type:	CW/pulsed 1470 nm
Power:	12W
Weight:	7.7lbs
Dimensions:	4" (h) x 9" (w) x 9" (d)

Proprietary Cooling Technology

Reduces unit size and weight considerably while providing similar power characteristics compared to other surgical lasers.

Energy delivered at 1470 nm

Highly absorbed by water and tissue allowing precise soft tissue treatment through its customized fiber optics.

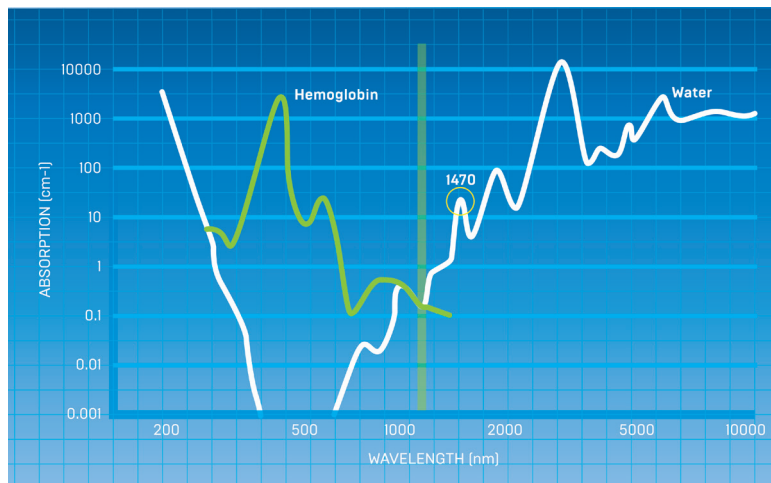

Specially designed fiber tip with small outer diameter

Short fiber tip for enhanced access into smaller and curved tracts, providing circumferential and uniform energy distribution.

- Forward firing tip with 360 degree array of energy
- One fiber per case
- Rounded tip for safety

Dedicated laser marking

Ensures precise positioning of the probe within the tract.



FORTEC medical

>> Start FORTEC Mobile Service Today

www.fortecmedical.com | 800.963.7101